

AD

TECHNICAL REPORT ARCCB-TR-94036

**EXPERIMENTAL METHODS IN RESIDUAL STRESS  
MEASUREMENT USING A POSITION-SENSITIVE  
SINGLE-EXPOSURE SCINTILLATION  
DETECTION SYSTEM**

**S.L. LEE  
M. DOXBECK  
G.P. CAPSIMALIS**

**DTIC  
ELECTE  
JAN 13 1995  
S G D**

SEPTEMBER 1994



**US ARMY ARMAMENT RESEARCH,  
DEVELOPMENT AND ENGINEERING CENTER  
CLOSE COMBAT ARMAMENTS CENTER  
BENÉT LABORATORIES  
WATERVLIET, N.Y. 12189-4050**



**APPROVED FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED**

**DTIC QUALITY INSPECTION**

**19950112 059**

#### DISCLAIMER

The findings in this report are not to be construed as an official Department of the Army position unless so designated by other authorized documents.

The use of trade name(s) and/or manufacturer(s) does not constitute an official indorsement or approval.

#### DESTRUCTION NOTICE

For classified documents, follow the procedures in DoD 5200.22-M, Industrial Security Manual, Section II-19 or DoD 5200.1-R, Information Security Program Regulation, Chapter IX.

For unclassified, limited documents, destroy by any method that will prevent disclosure of contents or reconstruction of the document.

For unclassified, unlimited documents, destroy when the report is no longer needed. Do not return it to the originator.

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE September 1994		3. REPORT TYPE AND DATES COVERED Final
4. TITLE AND SUBTITLE EXPERIMENTAL METHODS IN RESIDUAL STRESS MEASUREMENT USING A POSITION-SENSITIVE SINGLE-EXPOSURE SCINTILLATION DETECTION SYSTEM			5. FUNDING NUMBERS AMCMS: 6111.02.H611.1 PRON: 1A12Z1CANMBJ	
6. AUTHOR(S) S.L. Lee, M. Doxbeck, and G.P. Capsimalis				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) U.S. Army ARDEC Benét Laboratories, SMCAR-CCB-TL Watervliet, NY 12189-4050			8. PERFORMING ORGANIZATION REPORT NUMBER ARCCB-TR-94036	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES) U.S. Army ARDEC Close Combat Armaments Center Picatinny Arsenal, NJ 07806-5000			10. SPONSORING / MONITORING AGENCY REPORT NUMBER	
11. SUPPLEMENTARY NOTES				
12a. DISTRIBUTION / AVAILABILITY STATEMENT  Approved for public release; distribution unlimited			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words) In the Denver D-1000-A single-exposure position-sensitive scintillation detection (PSSD) system, surface residual stress measurement in any chosen direction is based on crystalline plane spacings determined in two directions normal to the surface made simultaneously in a single psi tilt. This technical allows fast, noncontacting, and nondestructive biaxial stress analysis. In this report, system performance is characterized by studying the noise, gain, and diffraction peak profiles as a function of diode array element. A four-point bend experiment was performed to determine the elastic constant of the 211 plane of body-centered cubic martensitic steel. Residual stress measurements were performed in several steel specimens and compared to measurements made on a similar system at Pennsylvania State University. Local software development allowed the single-exposure PSSD to run in a multiple-exposure mode for improved accuracy in biaxial stress analysis.				
14. SUBJECT TERMS Residual Stress, Biaxial Stress, X-Ray Diffraction, Position-Sensitive Scintillation Detector, Single-Exposure Technique			15. NUMBER OF PAGES 27	
			16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT UNCLASSIFIED	18. SECURITY CLASSIFICATION OF THIS PAGE UNCLASSIFIED	19. SECURITY CLASSIFICATION OF ABSTRACT UNCLASSIFIED	20. LIMITATION OF ABSTRACT UL	

## TABLE OF CONTENTS

ACKNOWLEDGEMENTS .....	iii
INTRODUCTION .....	1
INSTRUMENTATION .....	1
SINGLE-EXPOSURE X-RAY DIFFRACTION STRESS ANALYSIS .....	2
MULTIPLE-EXPOSURE X-RAY DIFFRACTION STRESS ANALYSIS .....	2
EXPERIMENTAL METHODS .....	2
Surface Preparation .....	2
Detector Gain .....	3
Detector Noise .....	3
Diffraction Profile .....	3
Specimen Displacement Error .....	4
Peak Location Determination .....	4
Cross Correlation Technique .....	5
FOUR-POINT BEND CALIBRATION AND ELASTIC CONSTANT DETERMINATION .....	5
STRESS ANALYSIS OF M106 SHELL SAMPLE .....	6
CONCLUSION .....	6
REFERENCES .....	8

### Tables

1. Residual Stresses in a Shell Specimen .....	6
--	---

### List of Illustrations

1. Schematic of the D-1000-A position-sensitive scintillation detection system .....	9
2. Single-exposure PSSD residual stress determination .....	10
3. Drawing of electropolishing apparatus .....	11
4. Left and right detector gain for various integration times .....	12
5. Left and right detector for various tube currents .....	13
6. Left and right detector gain signals indicating shift in peak locations due to electronic drift .....	14

7. Effect of noisy gain on noise- and gain-corrected diffraction profiles .....	15
8. Average left and right detector noise versus integration time .....	16
9. Right detector diffraction profiles for zero stress 400 mesh iron powder .....	17
10. Left detector diffraction profiles for zero stress 400 mesh iron powder .....	18
11. Right detector diffraction profiles for ferritic four-point bend specimen .....	19
12. Left detector diffraction profiles for ferritic four-point bend specimen .....	20
13. Four-point bend verification of x-ray stress versus strain gauge stress .....	21
14. M106 shell section inside surface axial stress right detector diffraction profiles .....	22
15. M106 shell section inside surface hoop stress left detector diffraction profiles .....	23
16. M106 shell section inside surface axial stress by $\sin^2\psi$ technique .....	24

Accession For	
NTIS CRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification .....	
By .....	
Distribution /	
Availability Codes	
Dist	Avail and/or Special
A-1	

## **ACKNOWLEDGEMENTS**

The authors would like to thank Professor C.O. Ruud and Mark Jacobs of Pennsylvania State University for running the shell specimen and Ellen Fogarty of Benét Laboratories for preparing the manuscript for publication.

## INTRODUCTION

Residual stresses, including macro and micro stresses, are caused by the nonuniform plastic flow of previous operations. Most manufacturing processes involving thermomechanical and chemical treatment of components generate residual stresses. In weapon systems, artificial processes are used to induce advantageous residual stress distribution for retarding fatigue crack growth and improving fatigue strength of the component.

In recent years, advances have been made in nondestructive techniques and instrumentation employing x-ray diffraction, ultrasonics, neutron diffraction, and Barkhausen noise techniques for residual stress determination. Among these, x-ray diffraction has been established as one of the most accurate and reliable methods (refs 1-3). The D-1000-A x-ray stress analyzer, manufactured by Denver X-Ray Instruments, is a single-exposure position-sensitive scintillation detection (PSSD) system. It allows fast, noncontacting, nondestructive biaxial surface stress analysis in metals, alloys, and other polycrystalline materials. This report summarizes the experimental methods employed during the calibration of this instrument. Utilizing the four-point bend apparatus with a bar specimen made of ASTM A723 gun steel, x-ray stresses and strain gauge results were compared, and the elastic constant for the 211 plane of martensitic steel was determined. Single-exposure stress analysis was made for several specimens. Local software development has been made to use the single-exposure PSSD for multiple-exposure ( $\sin^2\psi$ ) x-ray diffraction stress analysis.

## INSTRUMENTATION

The D-1000-A x-ray stress analyzer is a prototype instrument designed after a Ruud-Barrett PSSD (refs 4,5). The instrument utilizes a miniature x-ray tube that fits inside a cylinder with a 94-mm bore diameter. The operating principles for this detection system, shown in Figure 1, consist of converting incident x-rays into light through cadmium-zinc sulfide scintillation coating; transporting this light to a detector package through fiber optics bundles; amplifying the light intensity through an image intensifier; and converting it into electronic signals through a reticon photodiode array. The PSSD exploits x-ray sources of about 100 to 200 watts and provides a shorter data collection time interval of any x-ray diffraction detection system.

The system currently uses a chromium target tube, whose characteristic radiation does not fluoresce any major alloys in steel. Chromium K- $\alpha$  radiation reflects from the 211 plane of the body-centered-cubic (BCC) martensitic steel at a  $2\theta$  angle of 156.41 degrees. The reticon diode arrays, which are solid state scanners, consist of a row of silicon photodiodes, each with an associated storage capacitor upon which the photocurrent is integrated. The device has sensors on 25.4 micron centers with a sensing aperture 430 microns wide. There are 512 diodes in each of the two arrays. The array dimensions are length = 25.4 micron  $\times$  512 = 1.3 cm and width = 430 micron or 0.43 mm. The detector, 1.3 cm in length, spans an angle of 21.2 degrees for a focusing circle of radius  $R_o = 35$  mm. To obtain the required sensitivity, the silicon diode array is cooled to approximately -20°C by a copper heat sink connected to a water jacket. The cylinder is hermetically sealed and flushed with nitrogen so that water vapor does not condense on the cooled parts.

An IBM AT computer is connected to the system for data acquisition, control, storage, and analysis. Since the PSSD produces analog signals, data have to go through an analog-to-digital converter and associated electronics. The system software package uses a number of algorithms to correct electronic and mechanical hardware fluctuations, as well as x-ray focusing errors. Specific data fitting algorithms are used for background correction and diffraction peak determination.

## SINGLE-EXPOSURE X-RAY DIFFRACTION STRESS ANALYSIS

In x-ray diffraction, stress is determined indirectly by the strain it produces. Hooke's law relates stress and strain by the stiffness and compliance moduli. Stress and strain are second rank tensors, while the stiffness and compliance moduli are fourth rank tensors. From symmetry relationships in most crystal structures, the number of independent nonzero constants in these matrices is greatly reduced. Since stress perpendicular to a free surface is zero, our interest focuses on biaxial stress analysis.

Characteristic x-ray radiations scattered from the successive crystalline planes undergo constructive interference if certain geometric conditions satisfying Bragg's law are met. The interplanar lattice d-spacing acts as a strain gauge. Diffraction peak location measurements can be used to determine both residual and loading stresses. As shown in Figure 2, stress is determined from the shift in diffraction peaks located on two opposite sides of the diffraction cone in a single psi tilt:

$$\sigma = \frac{E(S_2 - S_1)}{(1 + \nu) \cdot 4R_o \sin^2 \theta \sin(2\beta)} \quad (1)$$

where E is Young's modulus,  $\nu$  is Poisson's ratio,  $(S_2 - S_1)$  represents the peak shift between the right and left diffraction peaks,  $R_o$  is the radius of the focusing circle,  $\beta$  is the angle the incident x-ray makes with the specimen normal, and  $\theta$  is the Bragg angle of the unstressed metal (ref 6).

## MULTIPLE-EXPOSURE X-RAY DIFFRACTION STRESS ANALYSIS

In single-exposure analysis using PSSD, a single psi tilt gives 12 points on the  $\sin^2\psi$  plot. Multiple-exposure x-ray diffraction stress analysis depends on the measurements from multiple psi angles of inclination. Using the PSSD, N number of psi tilts gives 2N number of points on the  $\sin^2\psi$  plot for stress determination by the least squares fit of a line. Multiple-exposure stress analysis using PSSD has advantages over both the single-exposure technique using PSSD and the conventional multiple-exposure technique using a single detector. Local software development has been made to use the PSSD for multiple-exposure biaxial stress analysis.

## EXPERIMENTAL METHODS

### Surface Preparation

X-ray penetration depends on the x-ray wavelength and the absorption coefficient of the specimen. For chromium x-rays on steel samples, penetration is about twenty microns. To make an accurate residual stress determination, removal of the specimen surface affected by cold working, e.g., grinding and polishing, corrosion, oxidation, and grease and polymer coatings, is necessary. This is done by machine and hand polishing, metallographic polishing, and then electropolishing of the specimen surface. We experimented with current density settings, concentrations of the polishing bath, and used both lead and stainless steel mesh as the cathode in the electropolishing procedure. Figure 3 shows the electropolishing setup used in our experiments.



## Detector Gain

When a constant x-ray flux is incident on the detector, pixel-to-pixel differences in total detection efficiency give rise to gain variations. The variations in the scintillation coating thickness, scintillation conversion efficiency, variable amplification of the image intensifier, nonuniform sensitivities of the silicon diode, fiber optics variations, and coupling factors can contribute to the nonuniformity of the gain signals.

To make pixel-to-pixel efficiency correction, gain signals were obtained by irradiating a glass surface with the x-ray beam at the same geometry the stress measurements were to be made. Since glass is an amorphous solid, it should have a uniform diffraction intensity distribution versus angle  $2\theta$  in the far back reflection range.

To investigate the uniformity and stability of the gain signal, a log of gain signal intensity versus diode array element was plotted at various integration times as shown in Figure 4 and at various tube currents as shown in Figure 5. Large pixel-to-pixel intensity variations were observed, especially in the right detector. As shown in the following sections of this report, diffraction profiles have to be corrected for nonuniform gain variations. Figure 6 shows that electronic drift in the detectors causes a gain shift. Two traces of gain signals versus diode array elements were recorded at subsequent times using the same integration time, tube current, and voltage, and they were plotted on the same graph. Figure 7 shows the effect of the noisy gain on the noise- and gain-corrected diffraction profiles. Gain stability and statistics are essential in accurate peak location determination. In this work, gain signals were obtained in 5- to 10-second integration times. It was continuously and carefully monitored for possible electronic gain drift.

## Detector Noise

Possible light leaks and the electronic dark current both in the image intensifier and diode arrays contribute to the detector noise, which is detected when no x-rays are applied. The dark current in the diode's arrays increases with the temperature of the arrays. Detector noise signals were collected for both the left and right detectors in the PSSD and were used to correct for the gain and profile readings.

Figure 8 is a plot of the average noise variation with detector integration time for the left and right detectors in the PSSD. This plot indicates that the right detector has a noise offset of 87, and the left detector has a noise offset of 25. To equalize the noise offsets, potentiometers within the instrument have to be adjusted. At an integration time near forty seconds, the curve begins to show an exponential behavior. At an integration time less than ten seconds, which was used throughout this work, a linear relationship between noise and the integration time is observed. A formula was derived from this plot and system software was modified to normalize integration time differences for the noise, gain, and profile

$$Noise(t_2) = Noise(t_1) + (t_2 - t_1) \cdot Slope \quad (2)$$

where  $t_1$  and  $t_2$  are two different integration times.

## Diffraction Profile

The intensity of the diffraction profile is proportional to the intensity of the characteristic radiation, which is in turn linearly dependent on the tube current and nonlinearly on the tube high voltage. In this work, a tube voltage of 80 kilovolts and a current of 6 milliamps were used. Integration time for the diffraction profiles ranged between 2 and 10 seconds. Enough sampling time should be allowed without saturating the charge capacitor.

Diffraction peak profiles are noise-corrected and divided by the gain efficiency to normalize the pixel-to-pixel gain sensitivities. System software was modified to perform this operation.

$$P(i) = (Profile(i) - Noise(i))/(Gain(i) - Noise(i)) \quad (3)$$

To obtain accurate stress measurement, sufficient integration time should be allowed for the profile and the gain signal. Diffraction peak profile determination was found to be difficult for noisy gain.

### Specimen Displacement Error

If the irradiated diffraction volume is displaced from the center of the diffractometer, there is a relative peak shift between  $\psi = 0$  and  $\psi = \psi$  due to specimen displacement error. The Denver x-ray stress analyzer has a  $\Omega$ -goniometer geometry, and the error in  $2\theta$  location is given by Eq. (4):

$$\delta(\Delta 2\theta) = (\Delta 2\theta)_{\psi=0} - (\Delta 2\theta)_{\psi=\psi} = \left( \frac{360}{\pi} \right) (\Delta X) \cdot \cos \theta \left[ \frac{1}{R_o} - \frac{\sin \theta}{R \sin(\theta + \psi)} \right] \quad (4)$$

where

$$R = R_o \frac{\cos(\psi + (90^\circ - \theta))}{\cos(\psi - (90^\circ - \theta))}$$

In Eq. (4),  $R_o$  is the diffractometer radius,  $\Delta X$  is the specimen displacement error,  $\theta$  is Bragg's angle, and  $\psi$  is the psi tilt. When the specimen is too high,  $\Delta X$  is positive. When the specimen is too low,  $\Delta X$  is negative. In the  $\sin^2\psi$  plot, high specimen position contributes to negative slope and compressive stress state. Low specimen position contributes to positive slope or tensile stress state. Precision specimen alignment by a micrometer dial-gauge indicator or laser-sensor indicator is suggested. ASTM recommends zero stress standard of loose powered iron. In this work, stress-free 400 mesh iron powder was used to obtain zero stress correction and to check for specimen alignment. Diffraction peak profiles for zero stress powder are given in Figures 9 and 10.

### Peak Location Determination

Before the diffraction peaks can be located, the raw diffraction peak intensity profiles must be processed by several correction algorithms as provided by the PSSD system software (ref 7):

1. Incoherent bremsstrahlung background radiation subtraction.
2. Lorentz polarization absorption intensity corrections.
3. Planer instead of curved fiber optics surfaces were used. To correct for the fiber optics alignment in the location of the center of the x-ray sensitive fiber optics, a semi-empirical parabolic regression method was used.
4. Iterative correction for the uncertainty in the detector circle to specimen distance  $R_o$ .

Li and Sun (ref 8) studied contributions to the background radiation in x-ray stress analysis following derivations from diffraction physics. They proposed a nonlinear equation for background correction dependent on  $2\theta$  and the  $\psi$  tilt, which agreed well with the experimental background curve. Furthermore, we propose that the cross correlation technique be used to determine the diffraction peak shift for weak and noisy diffraction profiles that might occur in certain crystalline structures. In this work, a straight line background fitting algorithm was used. The user determines a range for the starting and ending pixel locations to make the background correction. The user can select the percentage of the upper peak intensity profile in the analysis for the parabolic fit. We found that 70 to 80 percent of the maximum peak intensity is good for the parabolic peak determination.

### Cross Correlation Technique

The cross correlation technique can be applied to determine diffraction peak shift, especially for weak and noisy diffraction profiles that might occur in certain crystalline structures. By cross correlating diffraction peak intensity versus diode array element signal at  $\psi = 0$  with signals at  $\psi = \psi_1$  and  $\psi = \psi_2$ , respectively, residual stress can be calculated from the diffraction peak shift. An alternate method is to cross correlate the left detector signal ( $I^{\psi_1}$ ) and the right detector signal ( $I^{\psi_2}$ ). The cross correlation function for each  $\psi$  tilt is given by

$$F_j^{\psi} = \sum_{i=1}^{512} I_i^{\psi_1} + I_{i-j}^{\psi_2} \quad (5)$$

where  $i, j$  are diode array element indices. From the location of the maximum in the cross correlation function, peak shift ( $S_2 - S_1$ ) in Eq. (1) can be obtained.

### **FOUR-POINT BEND CALIBRATION AND ELASTIC CONSTANT DETERMINATION**

A four-point bend instrument was used to compare the x-ray and strain gauge stress measurements. A strain gauge was attached to a 6 by 1 inch ferritic steel specimen made from ASTM A723 gun steel, and both tensile and compressive stresses were applied. In Figures 11 and 12, the raw and corrected peak profiles for the ferritic steel specimen are shown for the right and left detectors, respectively. In Figure 13, the x-ray stress is plotted versus strain gauge stress. The elastic stress constant was determined from the slope of the curve as follows:

$$\sigma = \left( \frac{E_{hkl}}{1+\nu} \right) (S_2 - S_1) \left( \frac{0.0254}{4R_o \sin^2 \theta \sin(2\beta)} \right)$$

$$\frac{E_{hkl}}{1+\nu} = \left( \frac{\sigma}{S_2 - S_1} \right) \left( \frac{1}{K} \right) \text{ where } K = \frac{0.0254}{4R_o \sin^2 \theta \sin(2\beta)} \quad (6)$$

Our experiment indicated that the constant  $E_{hkl}/(1+\nu) = 23.9 \cdot 10^6$  psi. The quantity was determined to be  $23.5 \cdot 10^6$  psi for the same specimen using the stress analyzer at the Materials Research Laboratory, Pennsylvania State University during our visit.

## STRESS ANALYSIS OF M106 SHELL SAMPLE

In Figure 14, the raw and corrected right detector peak profiles for the inside surface axial stress of an M106 shell section are shown. The left detector inside surface hoop stress profiles of an M106 shell section are shown in Figure 15. In Figure 16, the multiple-exposure technique was used for the M106 shell inside axial stress determination. Psi tilts at 12, 20, 25, and 30 degrees give pairs of data points (1,5), (2,6), (3,7), and (4,8) on the graph.  $D$  is the interplanar spacing of the planes at nonzero psi and  $D_0$  is the interplanar spacing of an unstressed metal. The  $D_0$  value was obtained from the literature (ref 9). Linear regression analysis was used to obtain the slope for residual stress determination. In Table 1, the stress results measured at Benét Laboratories are given along with measurements from Pennsylvania State University. Errors quoted represent data dispersion in successive measurements, including statistical error, systematic error, and error due to goodness of fit in the multiple-exposure analysis. Residual stress accuracies are expected to be approximately  $\pm 10$  Ksi. The measurements are in good agreement--Benét's values were slightly higher than Pennsylvania State's.

Table 1. Residual Stresses in a Shell Specimen

Stress Orientation	Benét Labs*	Penn State**	Benét Labs***
	Stress (Ksi)	Stress (Ksi)	Stress (Ksi)
inside hoop	-42.1 $\pm$ 3.3	-39.1 $\pm$ 3.1	-43.5 $\pm$ 2.6
inside axial	-43.4 $\pm$ 3.3	34.6 $\pm$ 1.8	-40.3 $\pm$ 3.9
outside hoop	-39.9 $\pm$ 1.4	-35.0 $\pm$ 1.6	-46.3 $\pm$ 3.2
outside axial	-39.1 $\pm$ 3.2	-35.3 $\pm$ 1.7	-43.4 $\pm$ 2.6

\*Denver D-1000-A Single-Exposure Stress Analyzer (local software modification)

\*\*Pennsylvania State University Single-Exposure Stress Analyzer

\*\*\*Denver D-1000-A Multiple-Exposure Measurements (local software development)

## CONCLUSION

Experimental methods using Denver D-1000-A single-exposure position-sensitive scintillation detection system (PSSD) are given in this report. Characteristics of noise, gain, and diffraction profiles were studied as a function of system parameters such as integration time, number of iterations, and x-ray tube current and voltage. A four-point bend experiment was performed for an ASTM A723 steel specimen. Elastic constant of the 211 plane of BCC steel was determined. Local software development was made to allow the single-exposure PSSD to run in a multiple-exposure mode. The technique improves accuracy in the  $\sin^2\psi$  biaxial stress analysis using the PSSD. Residual stress analyses were made for several ordnance steel specimens.

Suggested software enhancements include cross correlation technique for peak shift determination; capability to determine peak position by profile fit to Gaussian, Lorentzian, and Pearson VII profile types; layer removal stress corrections; and triaxial stress analysis for complete stress tensor determination.

Suggested system hardware improvements include: (1) uniform and stable gain by recoating the scintillation detector; (2) stable current and voltage power supply and associated electronics; (3) precision specimen alignment by micrometer dial-gauge or laser-sensor indicators; (4) automated psi motion control for multiple-exposure stress analysis; and (5) automated radiation shutter. By the application of an extended or mobile photodiode array and x-ray tubes of other wavelengths, Denver D-1000-A PSSD can be expanded to study other metals, alloys, and composite materials.

## REFERENCES

1. M.E. Hilley, "Residual Stress Measurement by X-Ray Diffraction," *SAE J784a*, SAE Information Report, Society of Automotive Engineers, Inc., 1980.
2. B.D. Cullity, *Elements of X-Ray Diffraction*, 2nd Edition, Addison-Wesley Publishing Company, Reading, MA, 1978.
3. I.C. Noyan and J.B. Cohen, *Residual Stress, Measurement by Diffraction and Interpretation*, Springer-Verlag, New York, 1987.
4. C.O. Ruud, "Position-Sensitive Detector Improves X-Ray Powder Diffraction," *Industrial Research and Development*, January 1983, pp. 84-87.
5. C.O. Ruud, "X-Ray Analysis and Advances in Portable Field Instrumentation," *Journal of Metals*, Vol. 31, No. 6, June 1979, pp. 10-15.
6. J.T. Norton, "X-Ray Stress Measurement by the Single-Exposure Technique," *Advances in X-Ray Analysis*, Vol. 11, Plenum Press, NY, 1968, pp. 401-410.
7. *D-1000-A Stress Analyzer Instruction Manual*, Denver X-Ray Instruments, Inc., Version 2.0, May 1987.
8. J. Li and Y. Sun, "Background Analysis and Its Correction for X-Ray Stress Measurements," in: *Residual Stresses in Science and Technology*, (E. Macherauch and V. Hauk, eds.), International Conference on Residual Stresses, Vol. 1, 1987.
9. I.C. Noyan, "Determination of the Unstressed Lattice Parameter for Triaxial Residual Stress Determination by X-Rays," Office of Naval Research, Technical Report No. 14, July 1984.

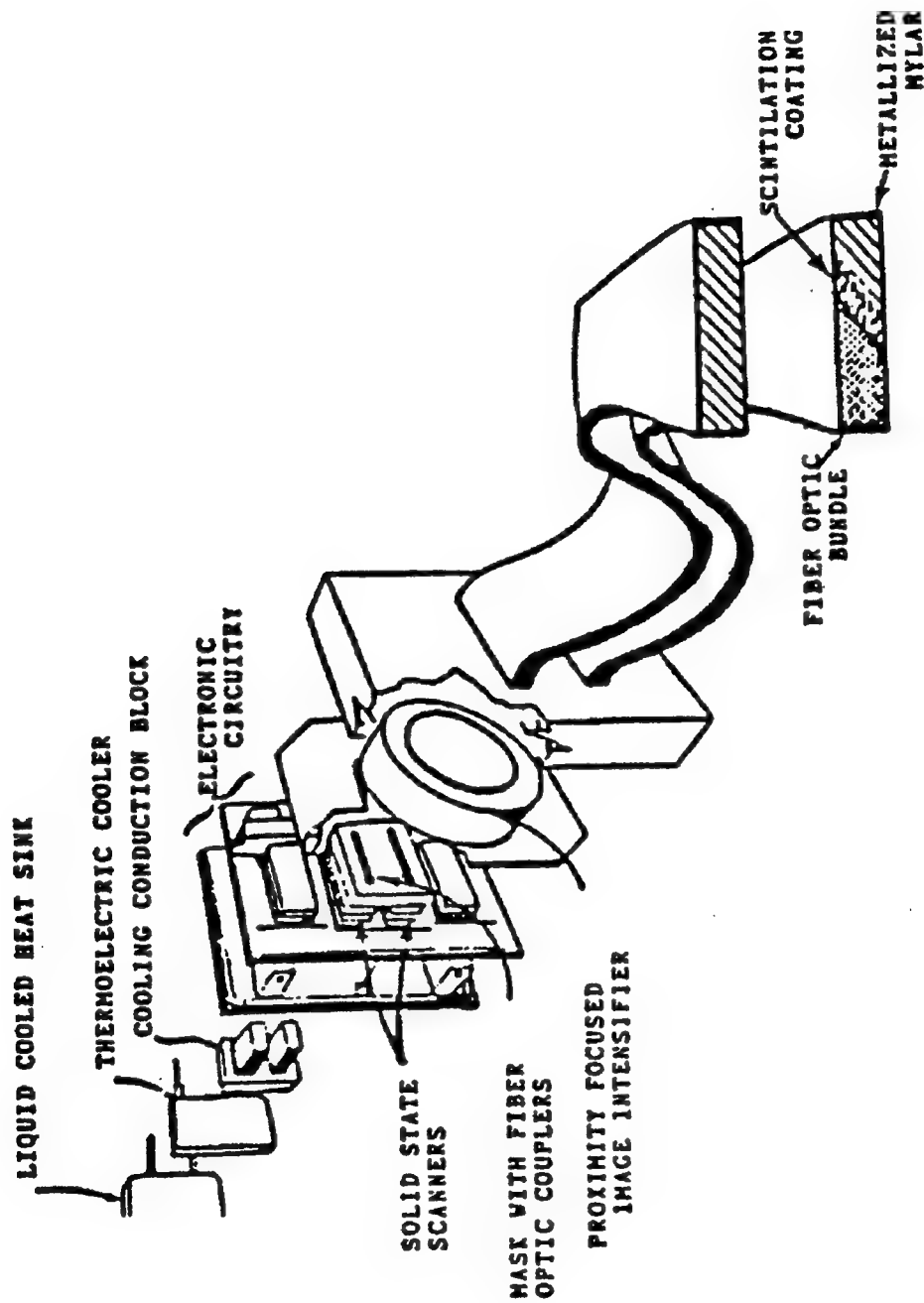


Figure 1. Schematic of the D-1000-A position sensitive scintillation detection system.





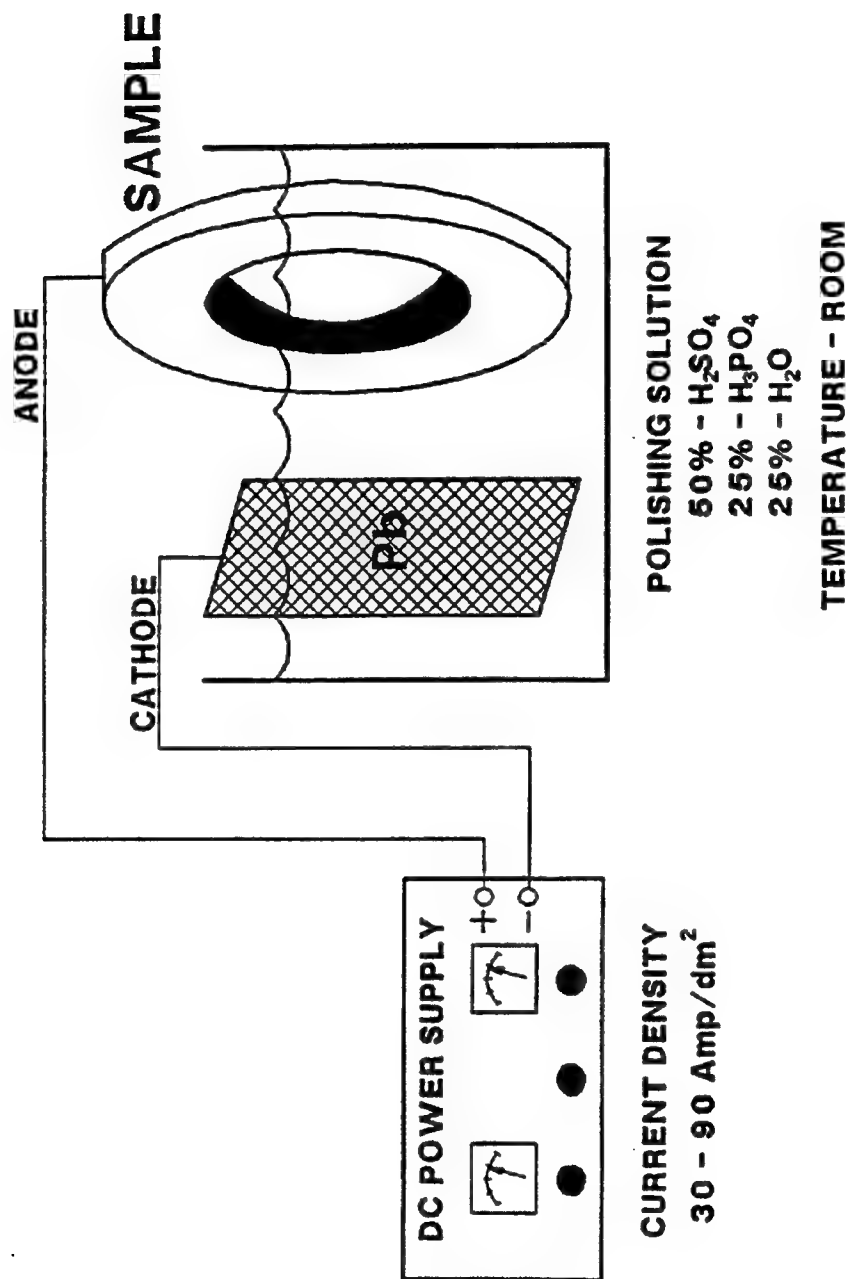


Figure 3. Drawing of electropolishing apparatus.

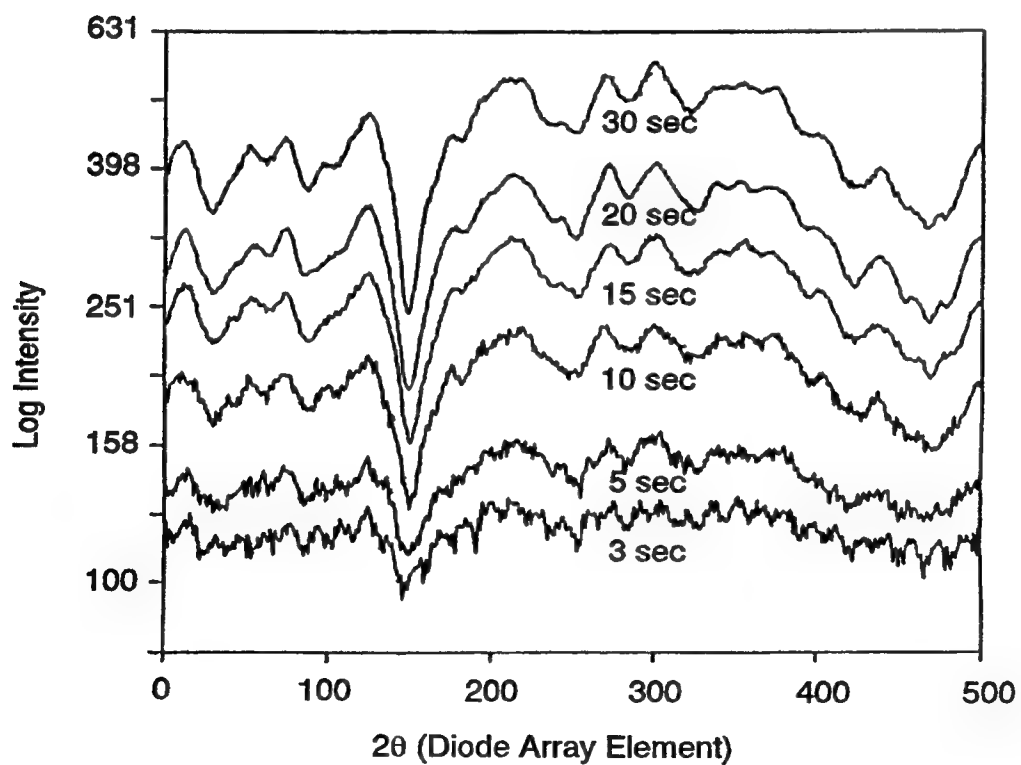
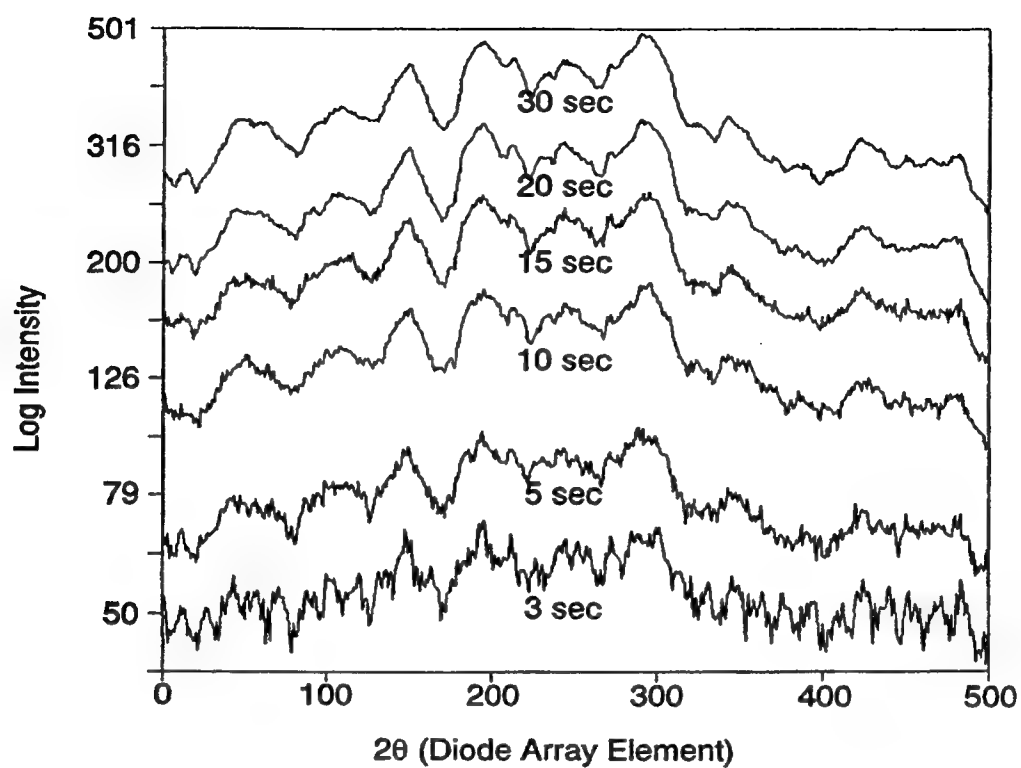


Figure 4. Left & right detector gain for various integration times.

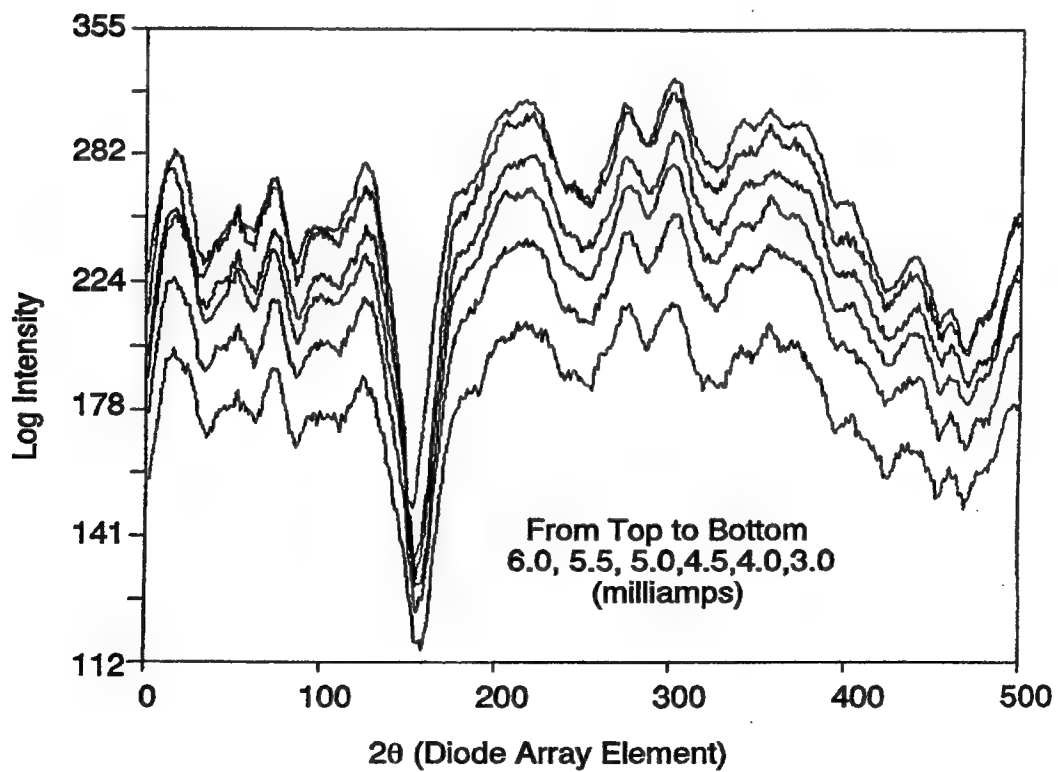
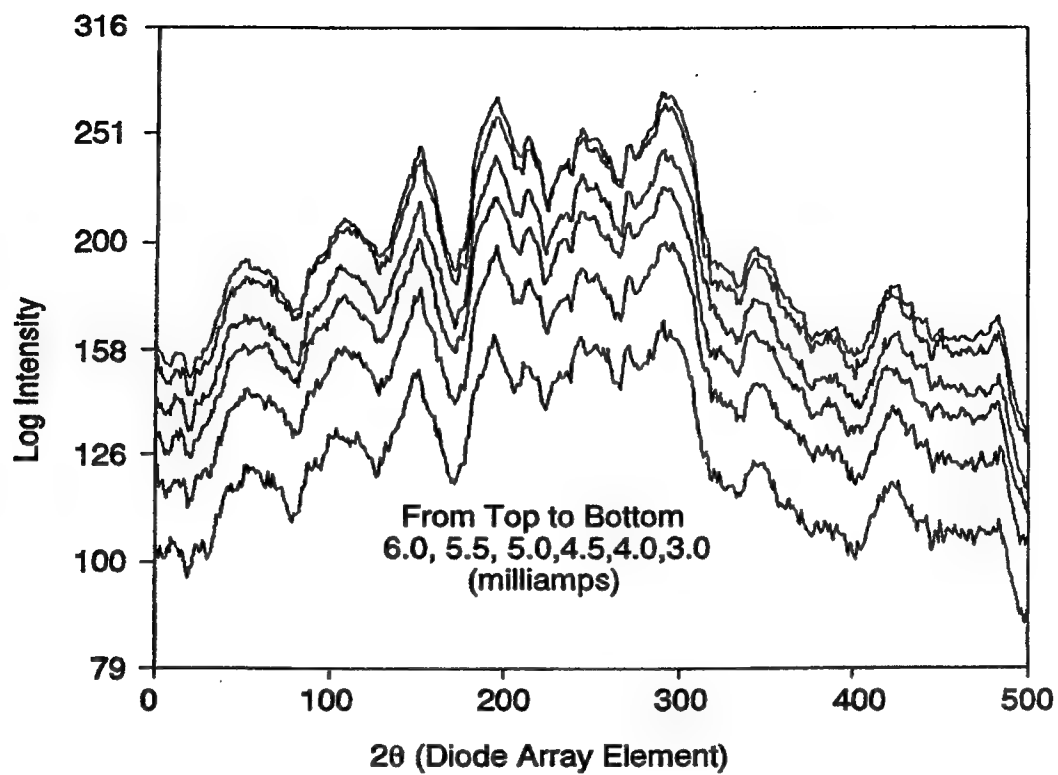


Figure 5. Left & right detector gain for various tube currents.

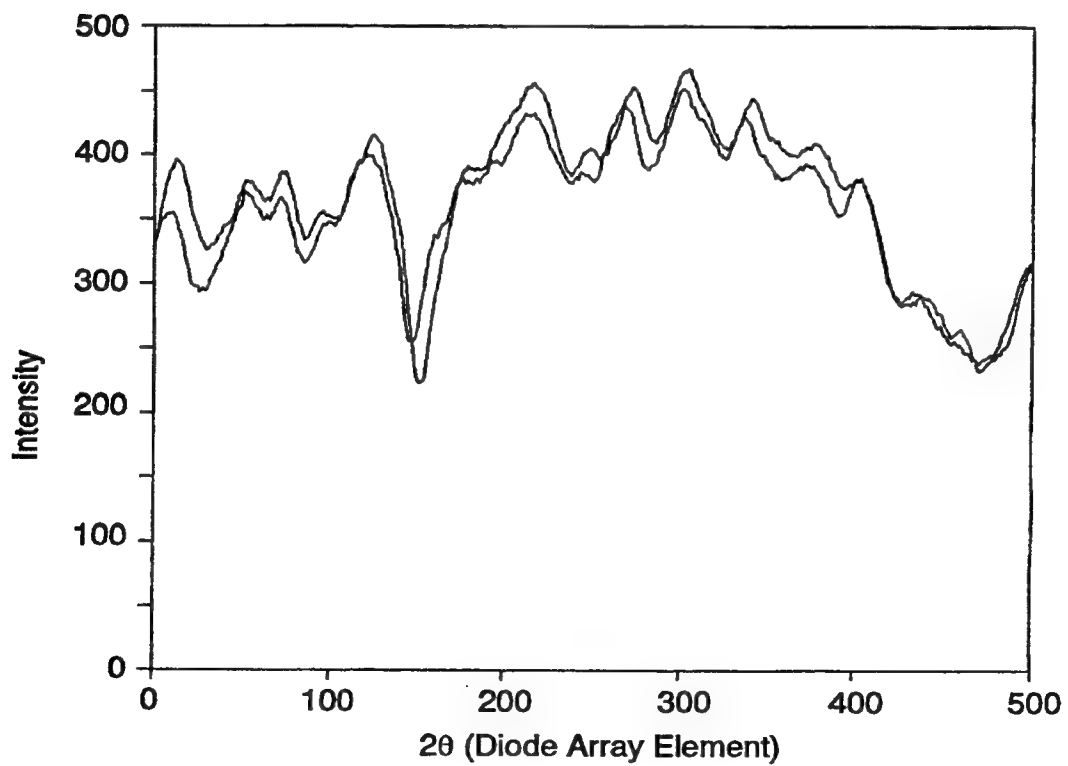
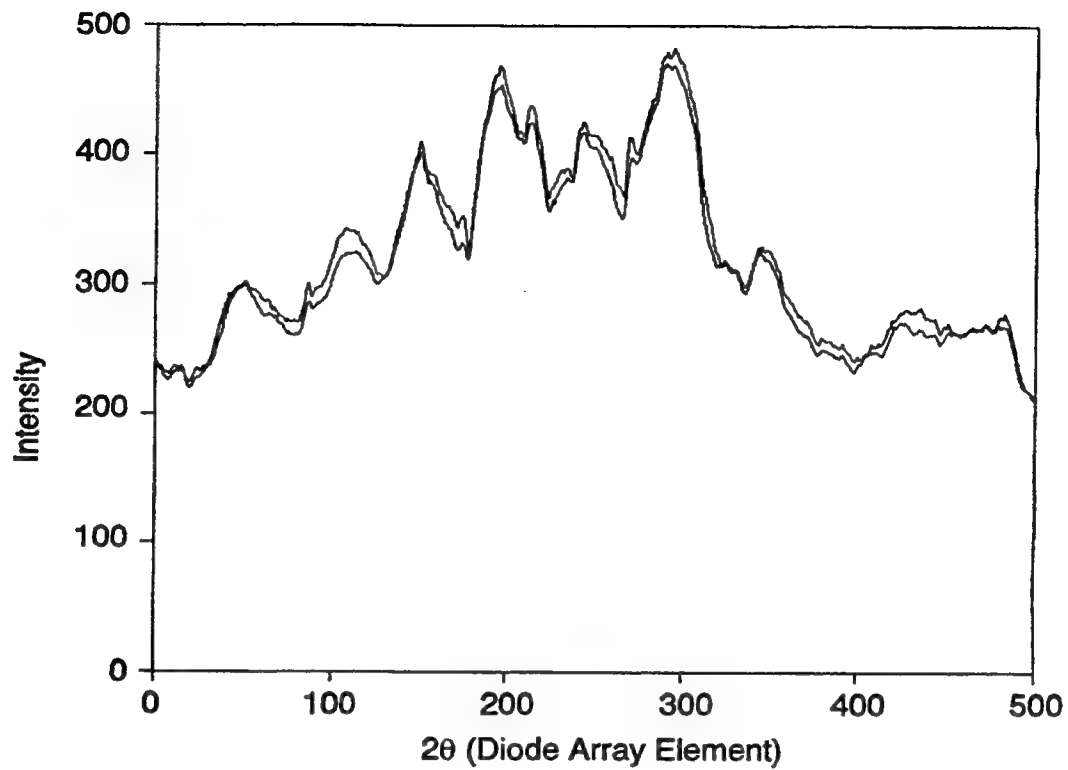


Figure 6. Left & right detector gain signals indicating shift in peak locations due to electronic drift.

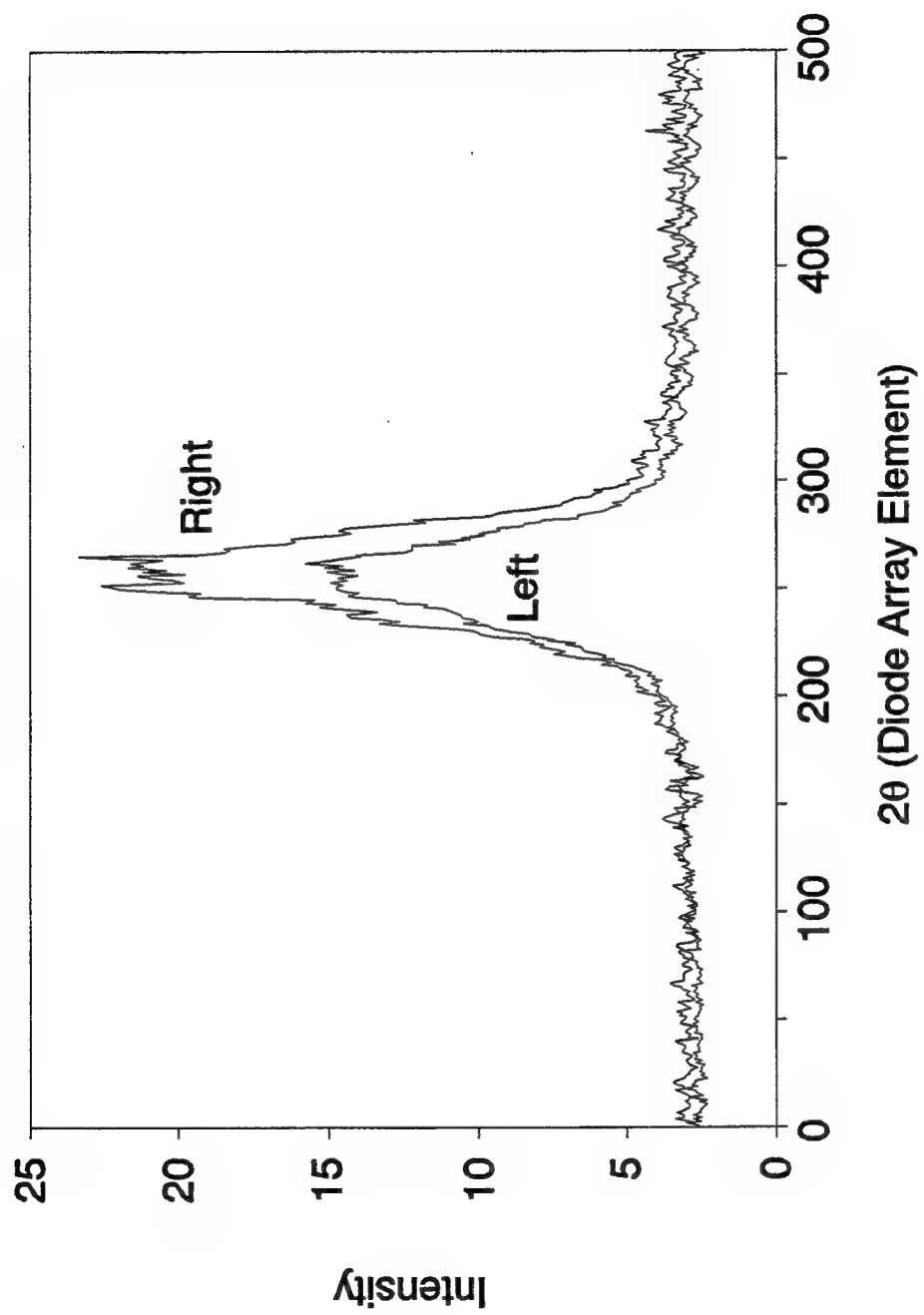


Figure 7. Effect of noisy gain on noise-and gain-corrected diffraction profiles.

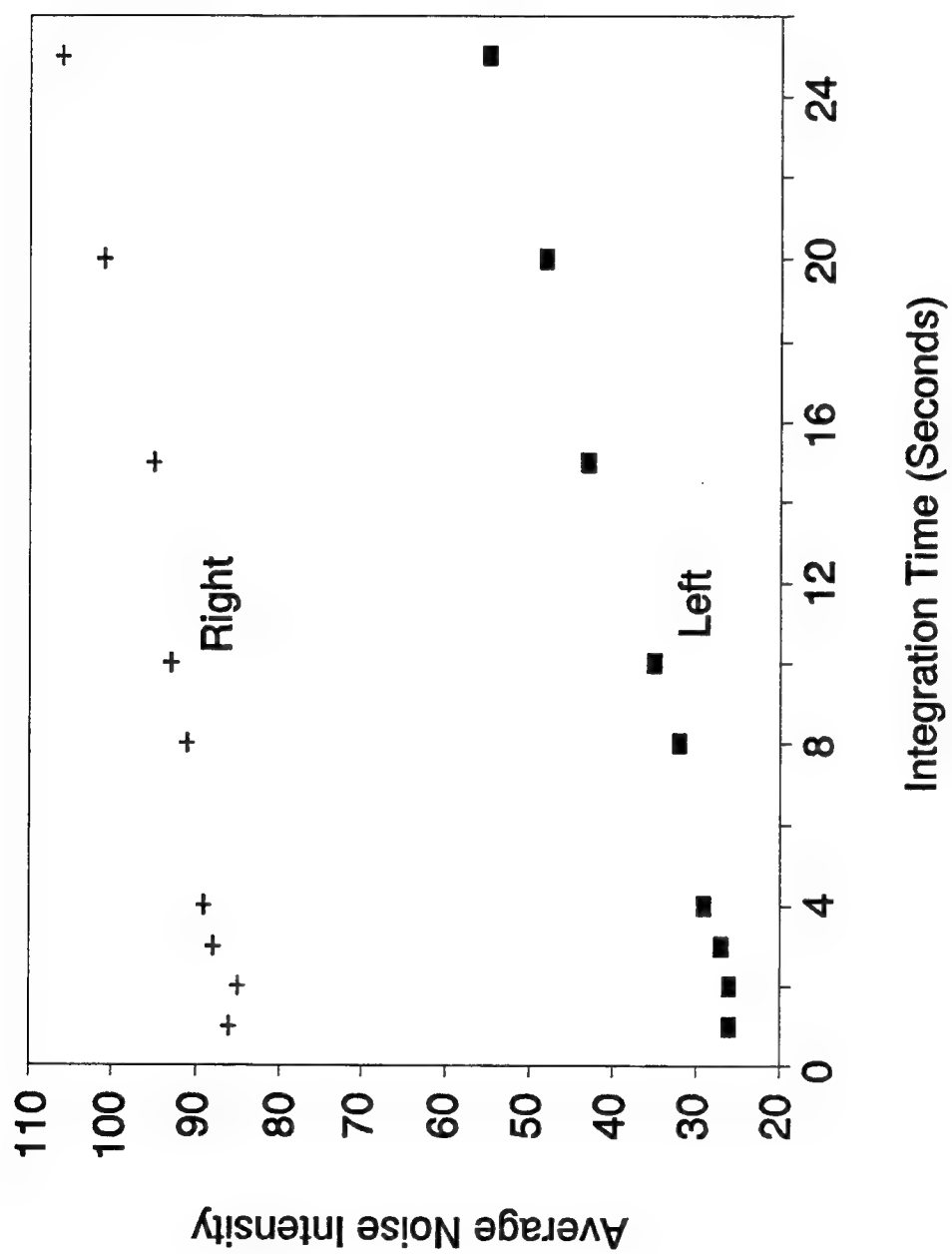


Figure 8. Average left and right detector noise versus integration time.

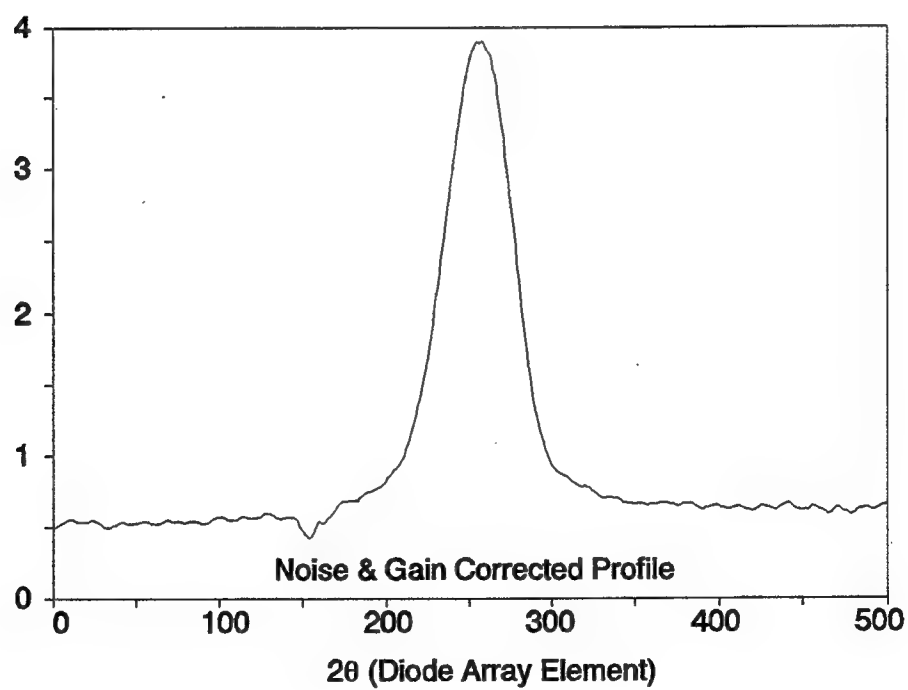
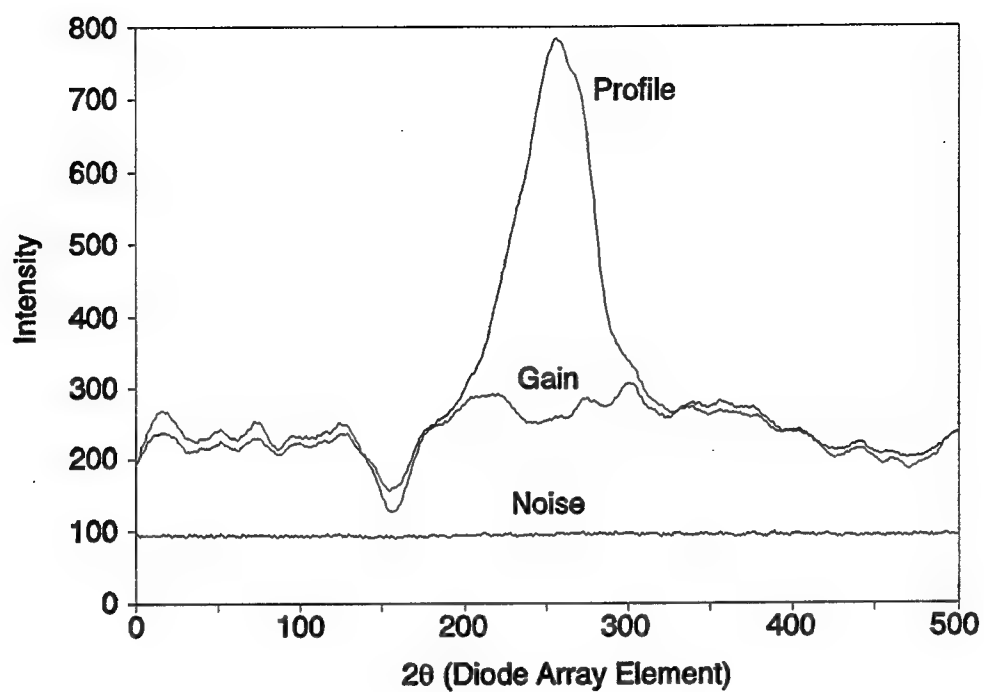


Figure 9. Right detector diffraction profiles for zero stress 400 mesh iron powder.

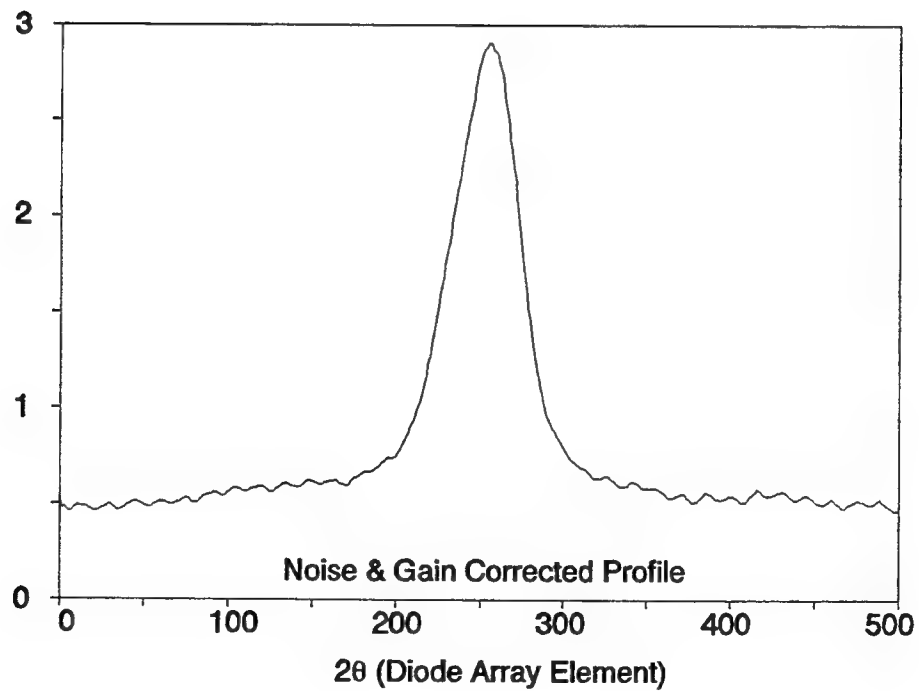
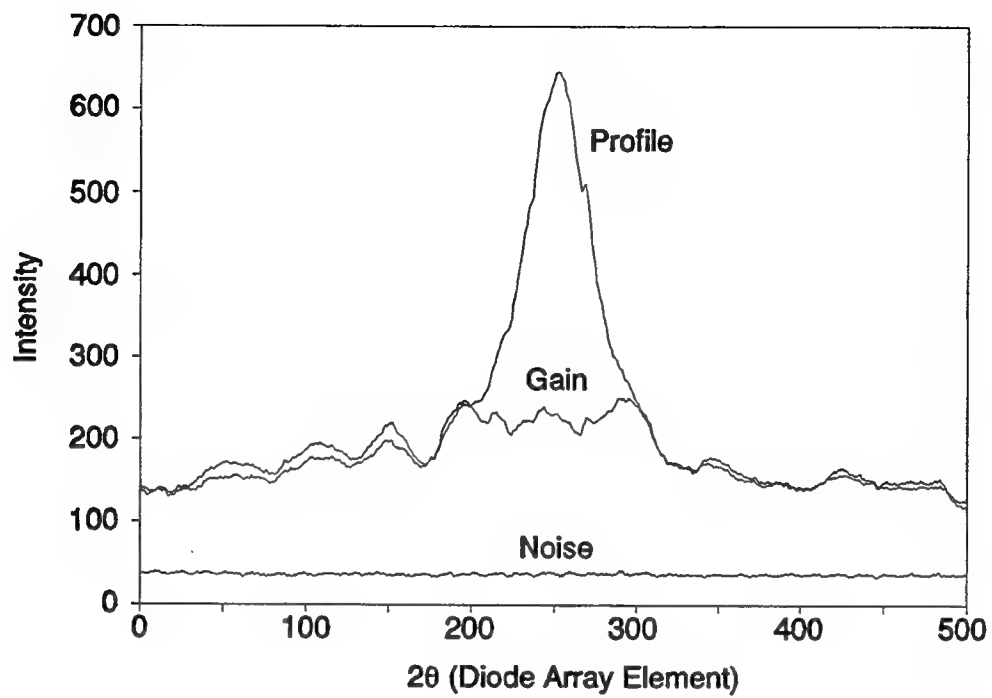


Figure 10. Left detector diffraction profiles for zero stress 400 mesh iron powder.



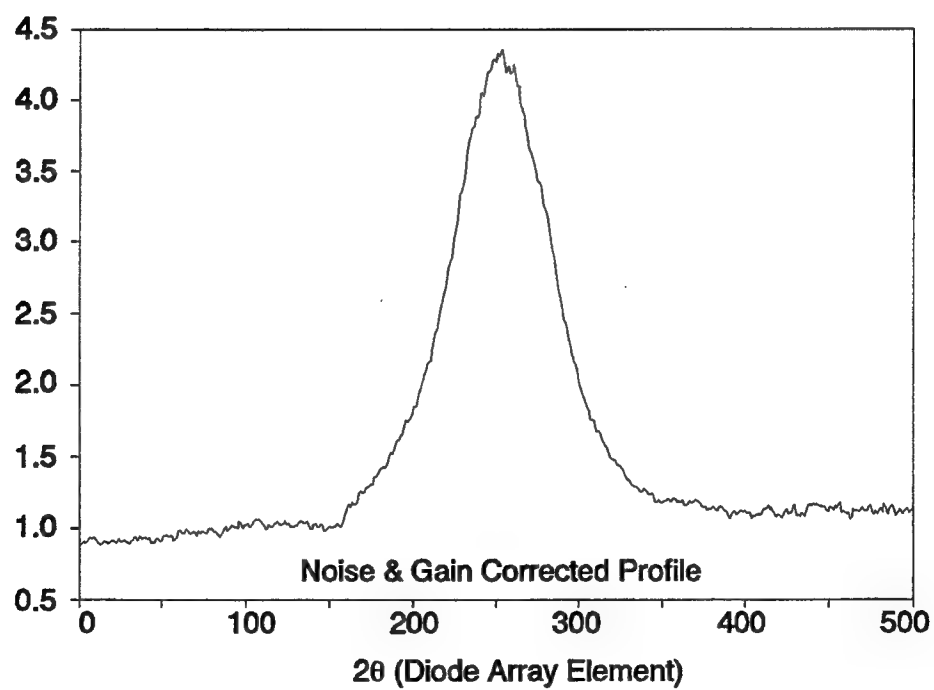
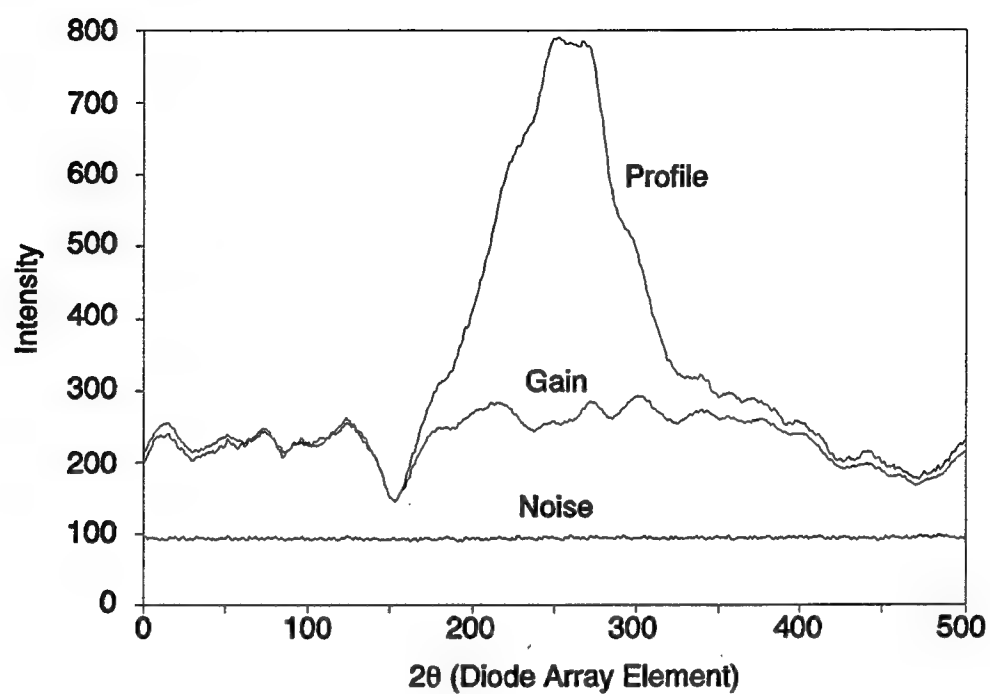


Figure 11. Right detector diffraction profiles for ferritic four-point bend specimen.

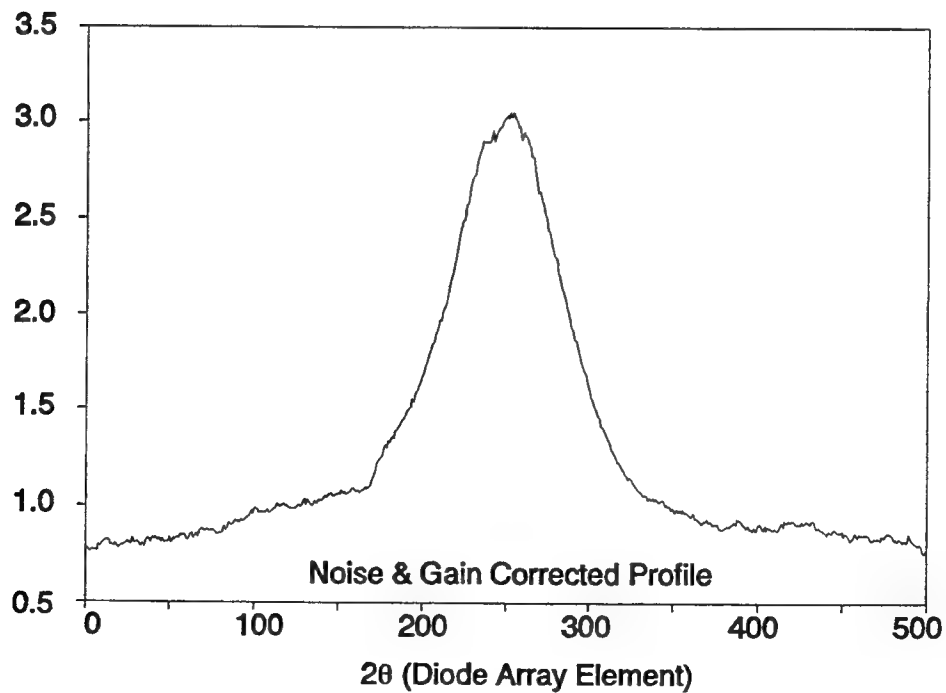
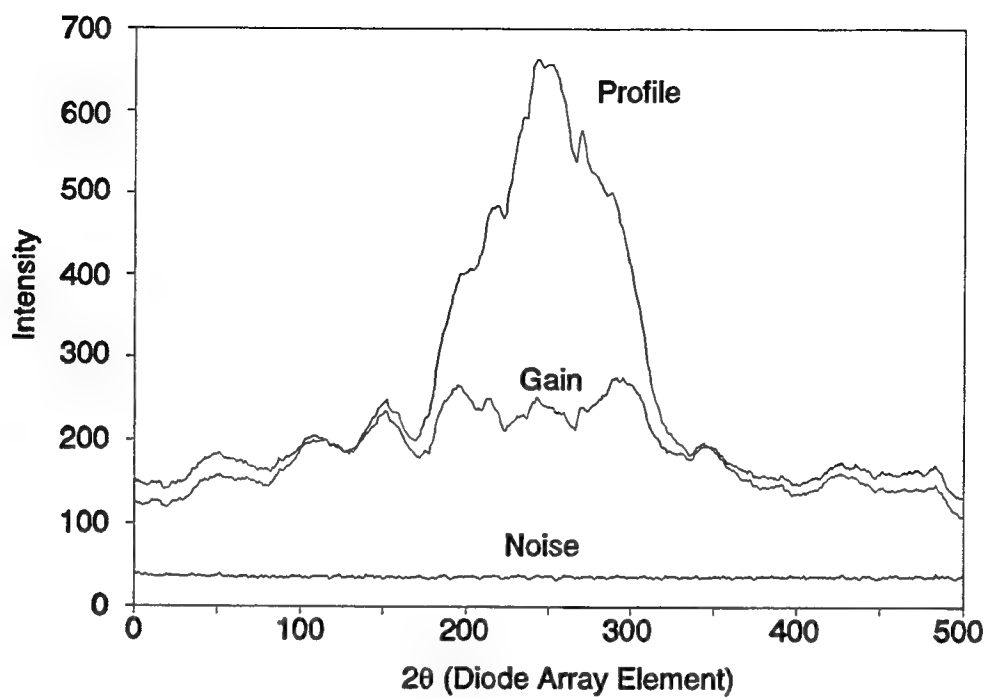


Figure 12. Left detector diffraction profiles for ferritic four-point bend specimen.

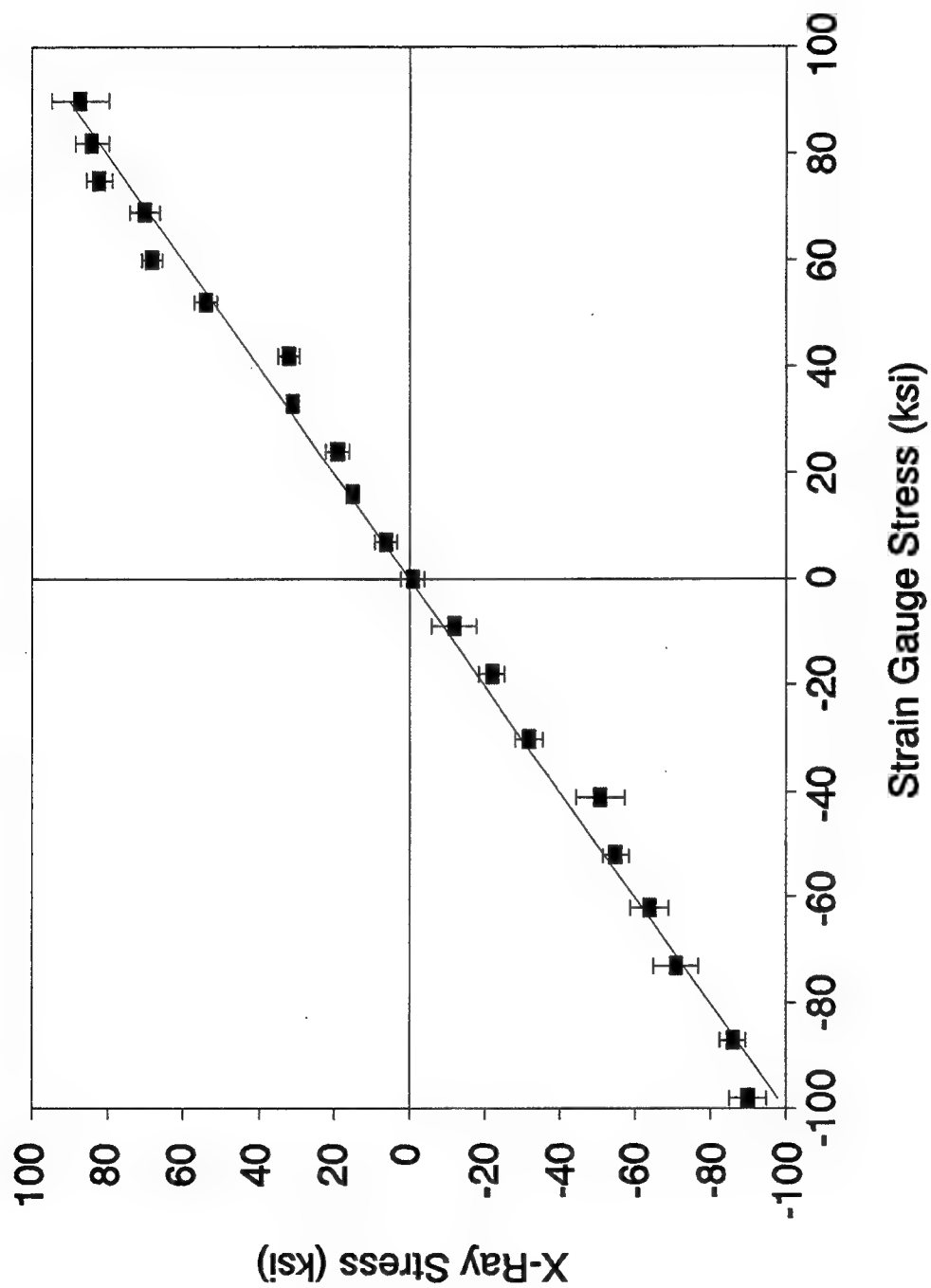


Figure 13. Four-point bend verification of x-ray stress versus strain gauge stress.

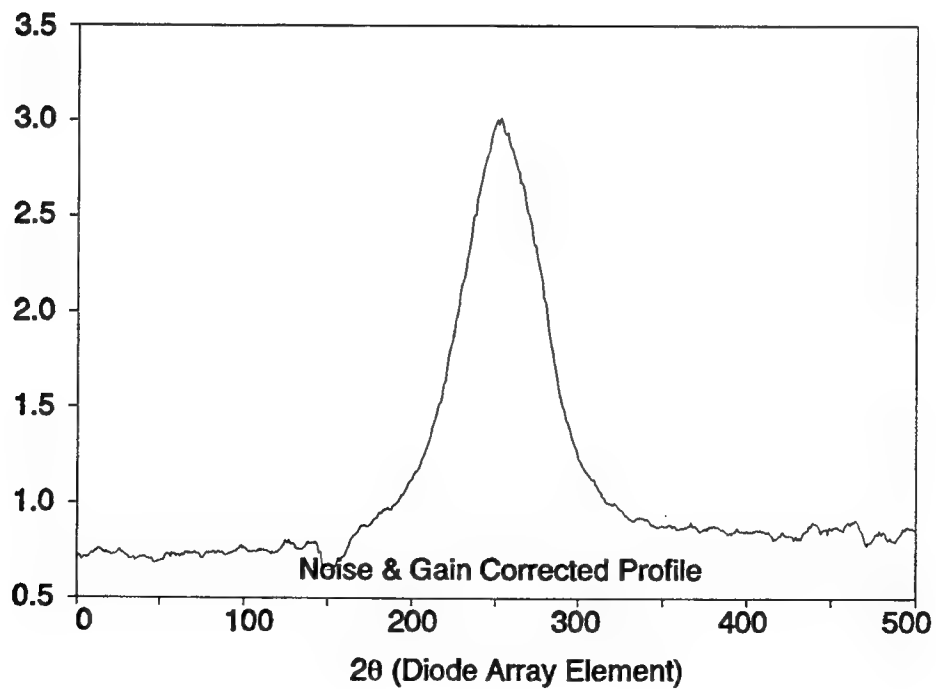
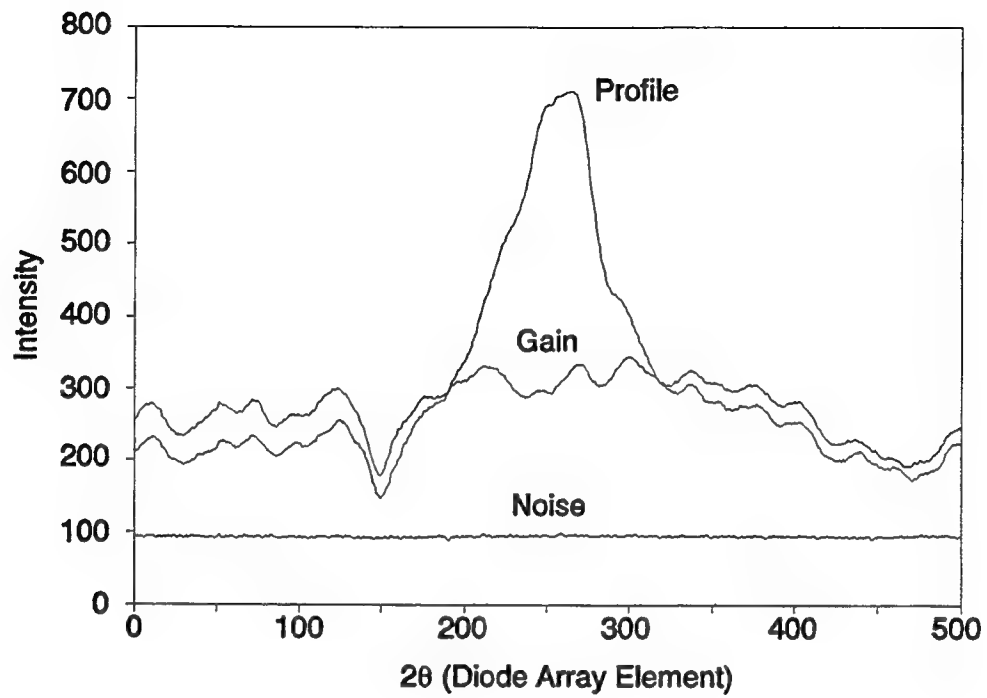


Figure 14. M106 shell section inside surface axial stress right detector diffraction profiles.

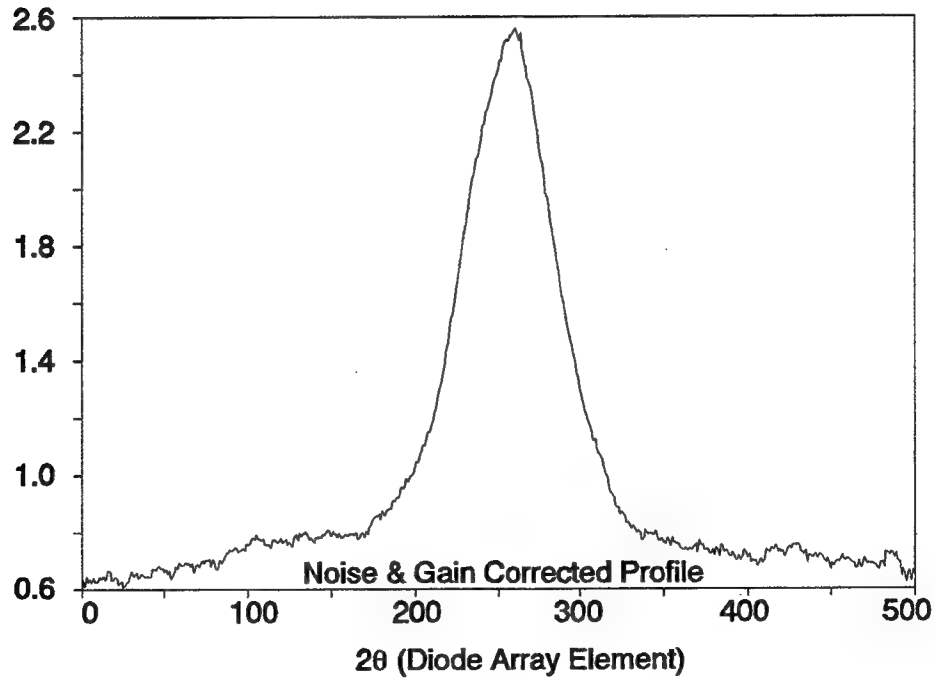
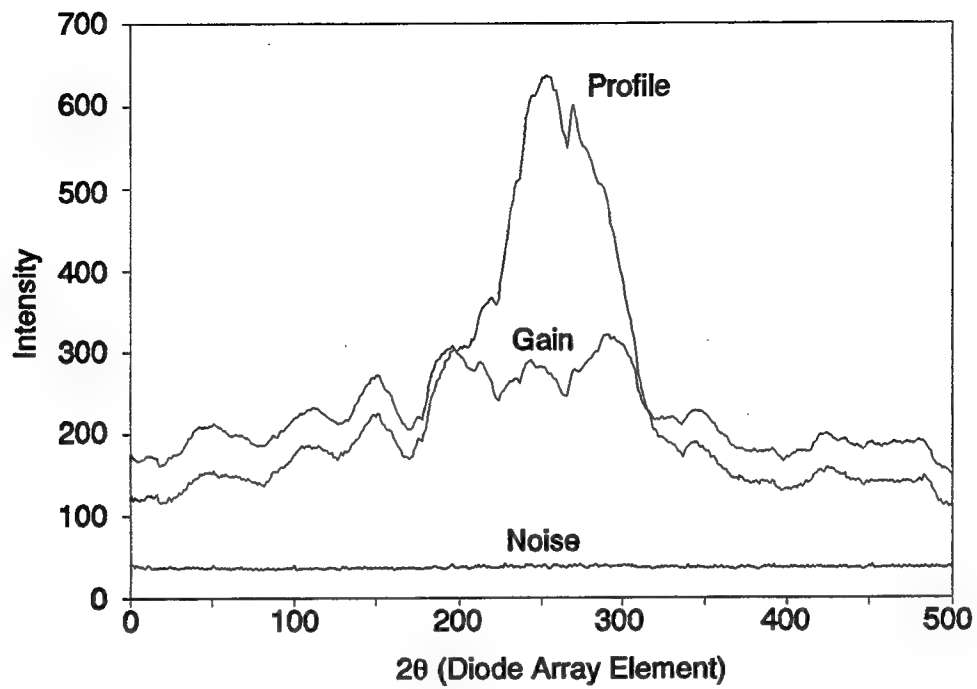


Figure 15. M106 shell section inside surface hoop stress left detector diffraction profiles.

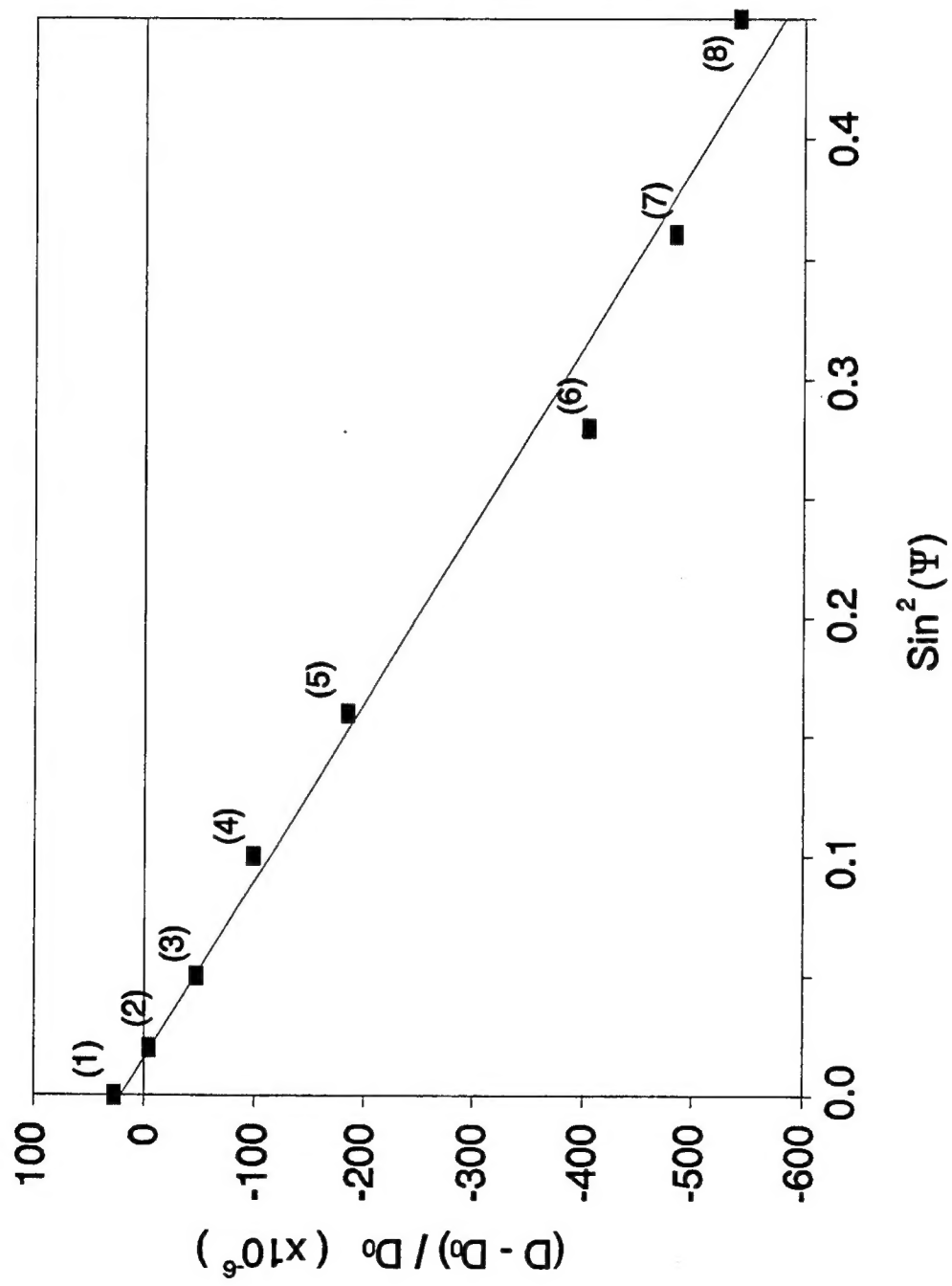


Figure 16. M106 shell section inside surface axial stress by  $\text{sin}^2(\Psi)$  technique.

---

TECHNICAL REPORT INTERNAL DISTRIBUTION LIST

	<u>NO. OF COPIES</u>
CHIEF, DEVELOPMENT ENGINEERING DIVISION	
ATTN: SMCAR-CCB-DA	1
-DC	1
-DI	1
-DR	1
-DS (SYSTEMS)	1
CHIEF, ENGINEERING DIVISION	
ATTN: SMCAR-CCB-S	1
-SD	1
-SE	1
CHIEF, RESEARCH DIVISION	
ATTN: SMCAR-CCB-R	2
-RA	1
-RE	1
-RM	1
-RP	1
-RT	1
TECHNICAL LIBRARY	
ATTN: SMCAR-CCB-TL	5
TECHNICAL PUBLICATIONS & EDITING SECTION	
ATTN: SMCAR-CCB-TL	3
OPERATIONS DIRECTORATE	
ATTN: SMCWV-ODP-P	1
DIRECTOR, PROCUREMENT & CONTRACTING DIRECTORATE	
ATTN: SMCWV-PP	1
DIRECTOR, PRODUCT ASSURANCE & TEST DIRECTORATE	
ATTN: SMCWV-QA	1

NOTE: PLEASE NOTIFY DIRECTOR, BENÉT LABORATORIES, ATTN: SMCAR-CCB-TL OF ADDRESS CHANGES.

---

---

TECHNICAL REPORT EXTERNAL DISTRIBUTION LIST

	<u>NO. OF COPIES</u>		<u>NO. OF COPIES</u>
ASST SEC OF THE ARMY RESEARCH AND DEVELOPMENT ATTN: DEPT FOR SCI AND TECH THE PENTAGON WASHINGTON, D.C. 20310-0103	1	COMMANDER ROCK ISLAND ARSENAL ATTN: SMCRI-ENM ROCK ISLAND, IL 61299-5000	1
ADMINISTRATOR DEFENSE TECHNICAL INFO CENTER ATTN: DTIC-FDAC CAMERON STATION ALEXANDRIA, VA 22304-6145	12	MIAC/CINDAS PURDUE UNIVERSITY P.O. BOX 2634 WEST LAFAYETTE, IN 47906	1
COMMANDER U.S. ARMY ARDEC ATTN: SMCAR-AEE	1	COMMANDER U.S. ARMY TANK-AUTMV R&D COMMAND ATTN: AMSTA-DDL (TECH LIBRARY) WARREN, MI 48397-5000	1
SMCAR-AES, BLDG. 321	1	COMMANDER U.S. MILITARY ACADEMY ATTN: DEPARTMENT OF MECHANICS	1
SMCAR-AET-O, BLDG. 351N	1	WEST POINT, NY 10966-1792	
SMCAR-FSA	1		
SMCAR-FSM-E	1		
SMCAR-FSS-D, BLDG. 94	1		
SMCAR-IMI-I, (STINFO) BLDG. 59	2	U.S. ARMY MISSILE COMMAND REDSTONE SCIENTIFIC INFO CENTER	2
PICATINNY ARSENAL, NJ 07806-5000		ATTN: DOCUMENTS SECTION, BLDG. 4484 REDSTONE ARSENAL, AL 35898-5241	
DIRECTOR U.S. ARMY RESEARCH LABORATORY ATTN: AMSRL-DD-T, BLDG. 305 ABERDEEN PROVING GROUND, MD 21005-5066	1	COMMANDER U.S. ARMY FOREIGN SCI & TECH CENTER ATTN: DRXST-SD 220 7TH STREET, N.E. CHARLOTTESVILLE, VA 22901	1
DIRECTOR U.S. ARMY RESEARCH LABORATORY ATTN: AMSRL-WT-PD (DR. B. BURNS) ABERDEEN PROVING GROUND, MD 21005-5066	1	COMMANDER U.S. ARMY LABCOM MATERIALS TECHNOLOGY LABORATORY ATTN: SLCMT-IML (TECH LIBRARY) WATERTOWN, MA 02172-0001	2
DIRECTOR U.S. MATERIEL SYSTEMS ANALYSIS ACTV ATTN: AMXSY-MP ABERDEEN PROVING GROUND, MD 21005-5071	1	COMMANDER U.S. ARMY LABCOM, ISA ATTN: SLCIS-IM-TL 2800 POWER MILL ROAD ADELPHI, MD 20783-1145	1

---

NOTE: PLEASE NOTIFY COMMANDER, ARMAMENT RESEARCH, DEVELOPMENT, AND ENGINEERING CENTER, U.S. ARMY AMCCOM, ATTN: BENET LABORATORIES, SMCAR-CCB-TL, WATERVLIET, NY 12189-4050 OF ADDRESS CHANGES.

---



---

TECHNICAL REPORT EXTERNAL DISTRIBUTION LIST (CONT'D)

	<u>NO. OF COPIES</u>		<u>NO. OF COPIES</u>
COMMANDER		COMMANDER	
U.S. ARMY RESEARCH OFFICE		AIR FORCE ARMAMENT LABORATORY	
ATTN: CHIEF, IPO	1	ATTN: AFATL/MN	1
P.O. BOX 12211		EGLIN AFB, FL 32542-5434	
RESEARCH TRIANGLE PARK, NC 27709-2211			
DIRECTOR		COMMANDER	
U.S. NAVAL RESEARCH LABORATORY		AIR FORCE ARMAMENT LABORATORY	
ATTN: MATERIALS SCI & TECH DIV	1	ATTN: AFATL/MNF	1
CODE 26-27 (DOC LIBRARY)	1	EGLIN AFB, FL 32542-5434	
WASHINGTON, D.C. 20375			

NOTE: PLEASE NOTIFY COMMANDER, ARMAMENT RESEARCH, DEVELOPMENT, AND ENGINEERING CENTER, U.S. ARMY AMCCOM, ATTN: BENET LABORATORIES, SMCAR-CCB-TL, WATERVLIET, NY 12189-4050 OF ADDRESS CHANGES.

---